


CGM Pediatric Electric Bed

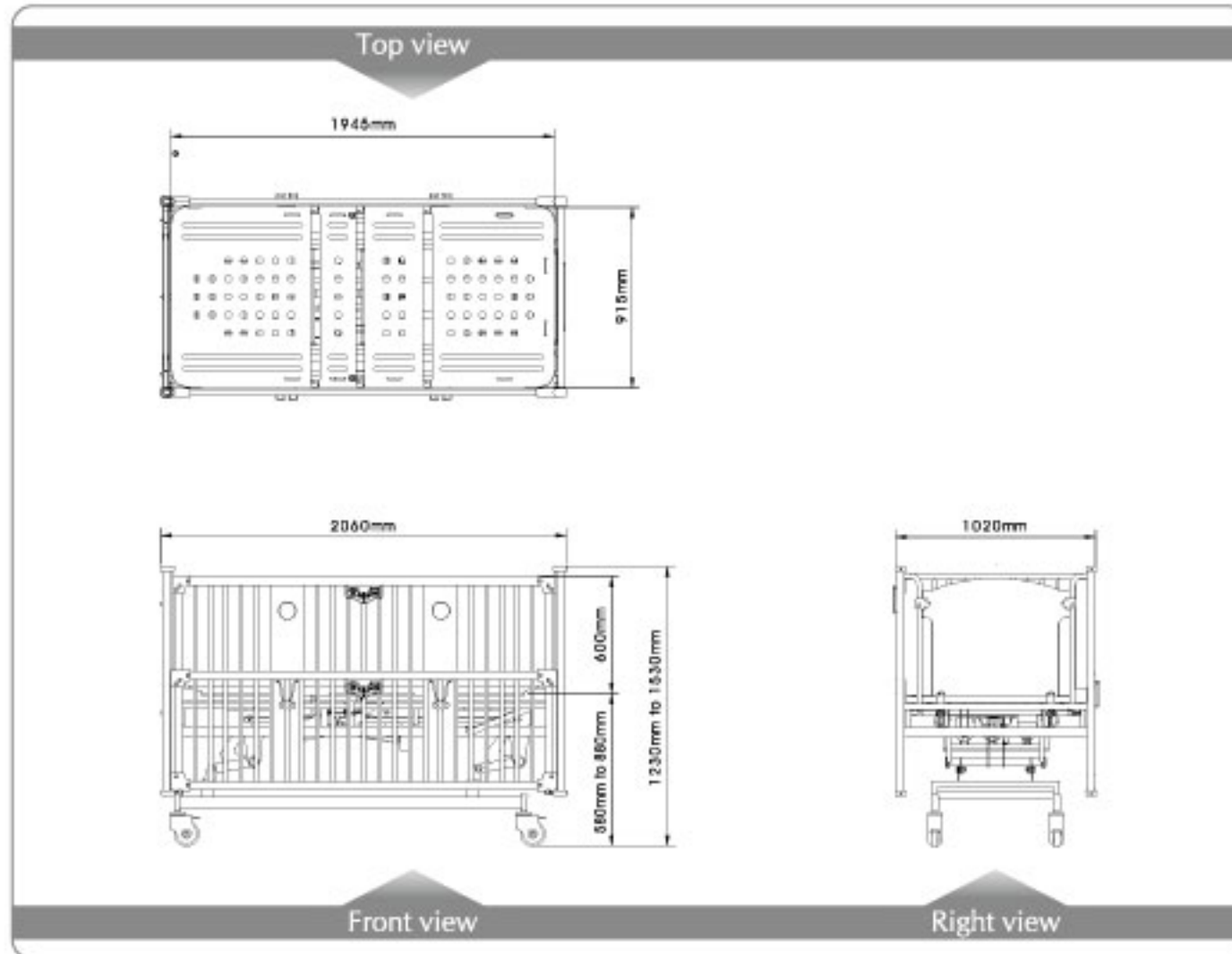

Specification

Dimension		Standard Accessories	
● Overall Size	L 206 x W 102 cm	● Fixed I.V. Pole	1PC
● Sleeping Surface	L 194 x W 91.5 cm	Optional Accessories	
● Siderail	L 196 x H 60 cm	● Mattress	
Positioning		● Oxygen Bottle Holder	
● Hi-Low Range	58 - 88 cm	● Under-Bed Basket	
● Backrest	0 - 75 °	● Charger Set	
● Knee Gatch	0 - 40 °	● Rechargeable Battery	
● Max. Weight Capacity	230 kg (500 lb)	● Weight Scale System	
Electric Standard		● Central Brake Caster	
● Voltage	220 ~ 240/100 ~ 120 VAC 50/60 Hz	● IV Pole	

CGMC reserves the right to make changes without notice in design and specification.



Chang Gung Medical Technology Co., Ltd.
Chang Gung Medical Supplies & Equipment Corp.

Office: 5F., No.201, Tung-Hwa North Rd., Taipei 105, Taiwan, R.O.C.
Tel: 886-2-87122948 Fax: 886-2-25140620

Factory: 2F., No.118, Nan Lin Rd., Taishan Shiang, Taipei 243, Taiwan, R.O.C.

E-mail: cgmc@cgmc.com.tw Http://www.cgmc.com.tw



CGEAP0908



Chang Gung Medical Technology Co., Ltd.
Chang Gung Medical Supplies & Equipment Corp.



Pediatric Electric Bed



1 Fixed I.V. Pole



- The stainless I.V. Pole can be placed on the head, the middle and the end of bed for hanging I.V. bottle and IVAC, according to medial care requirement.
- The height of the I.V. Pole is adjustable from 104 to 142 cm, measured from bed surface.

2 Siderail



- The fence-style siderail is total anti-rust and an unique design for complete security control.
- The safety height is to reach the standard of 60 cm to ensure kids' safety.

3 Parenting



A wide sleeping surface enables parents to sleep and accompany with their kids.

4 Caster

5* casters with dust-proof cover can prevent effectively dust from rolling inside the bearing.

5 Siderail Elevation Grip



- The smart grip is only for adults to move up/down and lock/unlock the siderail. By functional grips with 4-fixed position on both sides of tracks, the whole security system has been greatly upgraded.
- Gas-spring mechanism set on both sides of the siderails is to assist the siderail up and down easily.

6 Siderail ABS Cover

The ABS-made cover can keep the siderail from damages and increase its practicability. It also features warm color, easy-clean, smooth and comfortable touch.

7 Attendant Control Panel



- It is located at the foot-end of the bed.
- 3 rotation knobs for locking & unlocking the elevation of Backrest, Knee Gatch and Bed Hi-Low separately.
- 3 sets of leveling functions for adjusting the Backrest, Knee Gatch and Bed Hi-Low separately.

8 Embedded Control Panels



- 4 sets of Control Panels are embedded inside and outside of both siderails for elevation operations.
- 4 iconic indicators control the Backrest, Knee Gatch, Bed Hi-Low and Auto Movement elevations.

9 Electric Hi-Lo Functions



The Backrest, Knee Gatch and Bed Hi-Low are electric elevations, which can help the nursing to increase working efficiency and decrease occupation injury, ie. backache.

10 Removable ABS Platform

4-section air-permeable and washable ABS platform prevents bacteria from growing. The flat surface is very easy for cleaning and can prevent dirt accumulation.

11 Head Board

Pink wooden head board provides home-life impression. While emergency, the board can remove immediately.